

Medical Power Supply
Series GT*41134***, Series GT*96060*****,**
GTM96060-0706-1.0

Tested under

ANSI/AAMI ES 60601-1:2005 & A1:2012 & A2:2021 / CAN/CSA-C22.2 No. 60601-1:14 + A2:22 (R2022), Edition 3.2: Medical electrical equipment - Part 1: General requirements for basic safety and essential performance

CAN/CSA-C22.2 NO. 60601-1-6:11 + A1:15 + A2:21 (R2021) / IEC 60601-1-6:2020, Edition 3.2: Medical electrical equipment - Part 1-6: General requirements for basic safety and essential performance - Collateral standard: Usability

ANSI/AAMI HA 60601-1-11:15 + A1:2021 / CAN/CSA C22.2 NO. 60601-1-11:15 (R2020) + A1:21, Edition 2.1: Medical Electrical Equipment - Part 1-11: General requirements for basic safety and essential performance -- Collateral Standard: Requirements for medical electrical equipment and medical electrical equipment and medical electrical systems used in the home healthcare environment

File: E115461

MET Report: METC Class 137720

Approved: January 15, 2026

Applicant:

GlobTek, Inc.
186 Veterans Dr. Northvale,
NJ 07647,
USA

Prepared By:

Eurofins Electrical Testing Service (Shanghai) Co., Ltd.
Building 18, No. 2168 Chenhong Highway,
Minhang District, Shanghai,
China

For:

Eurofins Electrical and Electronic Testing North America, LLC
914 West Patapsco Avenue
Baltimore, Maryland 21230-3432
(410) 949-1802

This report shall not be reproduced except in full, without the express written consent of Eurofins Electrical and Electronic Testing North America, LLC

☒ MET Classification
☐ MET Recognition

☒ MET-C Classification
☐ MET-C Recognition

Table of Contents

CHANGE RECORD	3
DESCRIPTION	4
GENERAL REQUIREMENTS.....	8
MARKINGS.....	10
MANUAL/SERVICE INSTRUCTIONS	13
ALTERNATE LISTEE INFORMATION	14
APPLICANT'S RESPONSIBILITIES	15
PRODUCT MODIFICATIONS:	15
MANUFACTURING AND PRODUCTION-LINE TESTS AND DOCUMENTATION PERFORMED BY MANUFACTURER.....	17
CRITICAL COMPONENTS.....	18
CRITICAL DRAWINGS.....	25
FIGURES	26
ILLUSTRATIONS	36
TESTING CONSIDERATIONS	51
CONCLUSION	52

Change Record

Change Number	Description	Approval Date	Project Number	Amendment Engineer & Initial	Reviewer & Initials
--	--	--	--	--	--

Description

Product(s) Covered:

Medical power supply,

Model GTM96060-0706-1.0, Series GT*41134*****, GT*96060*****:

- The 1st “*” part can be ‘M’ or ‘-’ or ‘H’ for market identification and not related to safety.
- The 2nd “*” part can be “-” or “CC”, (“-” = Constant Voltage Model, CC = Constant Current Model)
- The 3rd “*” denotes the rated output wattage designation, which can be “01” to “06”, with interval of 1.
- The 4th “*” denotes the standard rated output voltage designation, which can be “03”, “03.3”, “04”, “06”, “12”, “15”, “18”, “24” or “36”.
- The 5th “*” is optional deviation, subtracted from standard output voltage, which can be “-0.1” to “-11.9” with interval of 0.1, or blank to indicate no voltage different.
- (The 4th “*” and 5th “*” together denote the output voltage, with a range of 3.3 - 36 volts.)
- The 6th “*” = Blank means directly plug in model series
 - = “-FW” means Class II open frame model with connector which is fixing on the PCB
 - = “-FWT2” means open frame model with appliance inlet with Class II inlet C8 respectively
- The last * denote any six character = 0-9 or A-Z or () [] or – or blank for marketing purposes.

Product Description:

- The equipment is AC-DC switching mode power supply for MED use.

Model Differences:

The different models are corresponding to four structure types respectively.

- First structure is direct plug-in power adapter with detachable plug portion, which is Class II apparatus. Two pieces of outer enclosure are enclosed with ultrasonic welding without screw.
- Second structure is open frame type. Interchangeable appliance inlets without earthing terminal can be mounted in medical device. The installation and use for the insulation construction shall be finally determined in the end product.
- Third structure Model GTM96060-0706-1.0, with rated output 5.0V, 1.5A, 7.5W. This structure is identical to the fifth structure except for different output ports and an additional inductor of output wire.
- Fourth structure is used in model series GT*41134***** and GT*96060*****, it uses F1 fuse in primary circuit and a LED indicator (optional) used in secondary circuit. GT*96060***** is identify with GT*41134***** except for model name.
- GT*41134*****, GT*96060*****, GTM96060-0706-1.0 direct plug-in series have same PCB; GT*41134*****-FW*, GT*41134*****-FWT2* open frame series are almost the same with direct plug-in series, but the PCB has 4 mounting holes.
- Transformers share the same construction, the turns may be added or reduced according different output voltage.

Description (Continued)

Electrical Rating:

100-240 Vac, Input for GT*41134****, GT*96060****, GTM96060-0706-1.0:

100-240V~, 50-60Hz, 0.3A or 0.6A

Output:

Model	Voltage	Max. current	Max. power
GT*41134**03*** GT*96060**03*** GT*41134**03.3*** GT*96060**03.3***	3.3V	1.8A	6W
GT*41134**04*** GT*96060**04***	3.4-4V	1.76A	6W
GT*41134**06*** GT*96060**06***	4.1-6V	1.46A	6W
GT*41134**12*** GT*96060**12***	6.1-12V	0.98A	6W
GT*41134**15*** GT*96060**15***	12.1-15V	0.50A	6W
GT*41134**18*** GT*96060**18***	15.1-18V	0.40A	6W
GT*41134**24*** GT*96060**24***	18.1-24V	0.33A	6W
GT*41134**36*** GT*96060**36***	24.1-36V	0.25A	6W
GTM96060-0706-1.0	5V	1.5A	7.5W

Description (Continued)

Engineering Considerations (Not For Field Representative's Use):

- The Medical power supplies have been investigated in accordance with the following standards:
ANSI/AAMI ES 60601-1:2005 & A1:2012 & A2:2021 / CAN/CSA-C22.2 No. 60601-1:14 + A2:22 (R2022), Edition 3.2: Medical electrical equipment - Part 1: General requirements for basic safety and essential performance;
CAN/CSA-C22.2 NO. 60601-1-6:11 + A1:15 + A2:21 (R2021) / IEC 60601-1-6:2020, Edition 3.2: Medical electrical equipment - Part 1-6: General requirements for basic safety and essential performance - Collateral standard: Usability;
ANSI/AAMI HA 60601-1-11:15 + A1:2021 / CAN/CSA C22.2 NO. 60601-1-11:15 (R2020) + A1:21, Edition 2.1: Medical Electrical Equipment - Part 1-11: General requirements for basic safety and essential performance -- Collateral Standard: Requirements for medical electrical equipment and medical electrical equipment and medical electrical systems used in the home healthcare environment.
- This product must be installed in accordance with all codes applicable to the location of the installation and in accordance with its instructions for use.
- Following essential performance(s), except those called out by particular standards, has been investigated and tested: None, which is documented in RMF N/A, along with test protocol N/A.
- Following applied part has been investigated in accordance with standards on the title page, including those identified by clause 4.6 of 60601-1 and documented in RMF: None.
- Following compartment of (or entire) the product has been investigated in accordance with standards on the title page for being used in Oxygen rich environment: None.
- Following supporting parts of the product has been investigated in accordance with standards on the title page: None with safety working load N/A kg and safety factor N/A.
- Following compartment of (or entire) the product has been investigated in accordance with standards on the title page for being used with flammable anesthetics: None.
- The insulation of plug was evaluated on compliance with two means of operator protection (2MOOP). The insulation system of the main equipment was evaluated on compliance with two means of patient protection (2MOPP).
- The product has been evaluated for a maximum ambient temperature of 50°C.
- The product is neither classified as Category AP nor APG.
- The product has been evaluated for use in a Pollution Degree 2 and overvoltage category II environment and a maximum altitude of 4000 m.
- The product is intended for indoor use only.
- For open frame model, temperature rise, leakage current and mechanical strength shall be reevaluated on this component when installed in the end product.
- The following clauses are not evaluated in this report:
Clause 11.7 Biocompatibility, referencing ISO 10993
Clause 17 EMC, referencing IEC 60601-1-2
- EUT is SMPS only, Basic safety is evaluated, no Essential Performance specified in manufacturer's RMF and the standards. Essential Performance is to be evaluated in end-product.

Description (Continued)

Note to Field Representative:

- A sample of each component listed below and a purchase order for the work described below at the current hourly rate shall be submitted to:*

Eurofins E&E NA, LLC
914 West Patapsco Avenue
Baltimore, Maryland 21230-3432

for reassessment processed under job #xxxxx for verification of construction against the associated drawings also listed below. The component(s) shall be subjected to an annual audit for continued compliance. The annual re-verification is a client incurred expense to be assessed at the current hourly rate at the time of the test. The estimated time for re-verification is also listed below.

Figure/ Item #	Component	Controlled Document Number	Re-Verification Type	Re-verification Maximum Estimated time (hours)
10-12/11	Transformer	Material Approval (Transformer)	Dimension - See Illustration 13-14, and Electric strength test	2

*Alternatively: If the evaluation is performed by the Eurofins E&E NA representative's lab other than the location above or by the MET representative during the Follow-up inspection, all data shall be returned to the Baltimore office listed above for surveillance tracking under the assigned job number mentioned above.

- The above inspections are a client incurred cost and will be billed at the hourly rate in place at the time of the inspection.

General Requirements

Scope of Requirements: The requirements contained within this section apply to all products contained within this Follow-Up Service Report File where applicable.

Definitions: (as defined or used in the context of the standard)

Term	Definitions
PCB:	Printed Circuit Board
Listed/Recognized Component:	A component evaluated to the applicable U.S. standards by a Nationally Recognized Testing Laboratory (NRTL).
Certified Component:	A component evaluated to the applicable Canadian standards by a Certification Organization (CO).
Listee:	Applicant

Measurements: All dimensions indicated in the body of this report are approximations unless otherwise indicated.

Corrosion Protection: All corrosive metals shall be provided with a means to protect from corrosion. Acceptable methods include painting, plating and galvanizing, or equivalent process. Dissimilar metals shall not be employed where reliable continuity is required.

Soldered Connections: All soldered connections shall be made mechanically secure before soldering. Tack soldering is not acceptable. Acceptable forms of mechanical securement include:

- A) Lead is inserted through an eyelet or opening of a terminal block prior to soldering.
- B) Lead is inserted into a U or V shaped slot in the terminal prior to soldering.
- C) Lead is wrapped around a terminal post prior to soldering.
- D) Lead is tied to adjacent lead with wire tie-wrap near termination point.

Electrical Connections: All electrical connections other than soldering shall be provided with positive detent, crimp type insulated Recognized Component connectors suitable for the voltage and temperatures involved. They shall be sized for the wire and mounting terminations. Where hazardous voltage or energy is involved, all wire connections to connectors shall employ a recognized method of double securement. Where fork-type lugs are used, they shall be snap-on or up-turned lug type.

Mechanical Assembly: All parts shall be secured by welding, bolts/nuts with lock or star washers, or thread forming screws.

Creepage and Clearances: See isolation diagram and table in illustration number x.

General Requirements (Continued)

Where present, the following items are required.


Internal Wiring: All internal wiring and connections are properly jacketed or enclosed within the equipment. Wiring is routed and secured to reduce the possibility of stress being transmitted to electrical connections, as necessary. All internal conductors in the secondary circuits are routed away from primary circuit conductors and from uninsulated live parts. There is no internal wiring subject to contact by the user when the product is employed as intended. The internal wiring is acceptable for conditions of service to which it will be subjected. Internal conductors consist of Recognized Component AWM insulated individual conductors; sized in accordance with the National Electric code and Canadian Electrical code, as may be applicable for the current expected in the conductor, rated VW-1, 300V, 105°C, and signal level ribbon wiring of flammability rating VW-1.

Interconnecting Cords and Cables: Flexible telecommunication cord and cable assemblies employed for interconnection between components are to be rated for and comply with temperatures, exposure to oil or grease and other conditions of service within the environment the product is to be utilized.

Markings

Etching, molding, die-stamping, silk-screening, stamped-, or etched-metal labels secured by rivets or screws are considered permanent. Recognized/Certified Component, Marking and Labeling Systems, and/or labels tested and deemed suitable for the surface to which it is applied is also considered permanent. Per the Canadian Electrical Code described in CSA C22.2 No. 0 General Requirements, Canadian product certification requires warning/cautionary markings in both English and French languages. It is the Applicant's responsibility to provide the listed Bilingual Markings shown below in accordance with the Canadian regulatory requirements. Each product is to be permanently marked with the following information:


- a. The MET Mark (refer to MET Applicant Contract), with the applicant/listee name or alternate listee name as identified within this report, trade name or trade mark, product model number, and a date of manufacture or serial number.
- b. Method of applying the MET Mark:
 - ☒ Direct Imprinting
 - ☐ Purchasing Labels from Eurofins E&E NA
 Approved MET Mark:

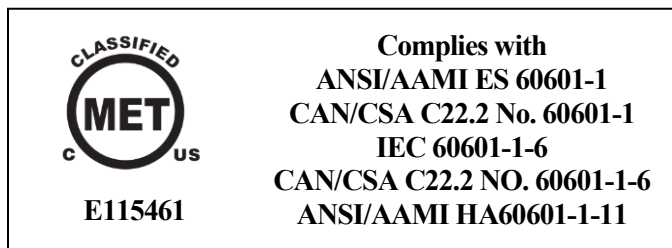


- c. A rating label adjacent to the inlet connector identifying the input voltage, current or power, frequency if applicable, replaceable fuse ratings if applicable (including voltage, current, interrupt capacity, rupturing speed) for the equipment.
- d. Near each terminal and control, there shall be a marking to identify the function.

Markings (Continued)


(Representative)

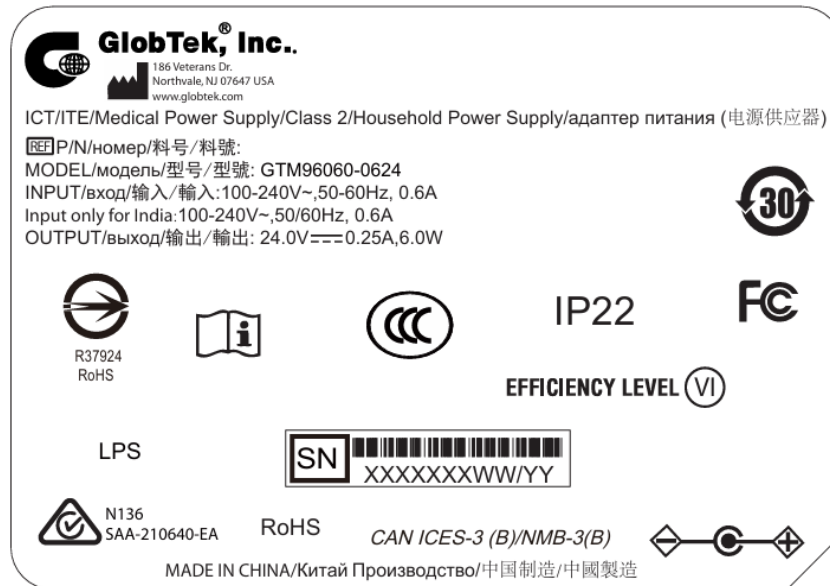

Markings (Continued)

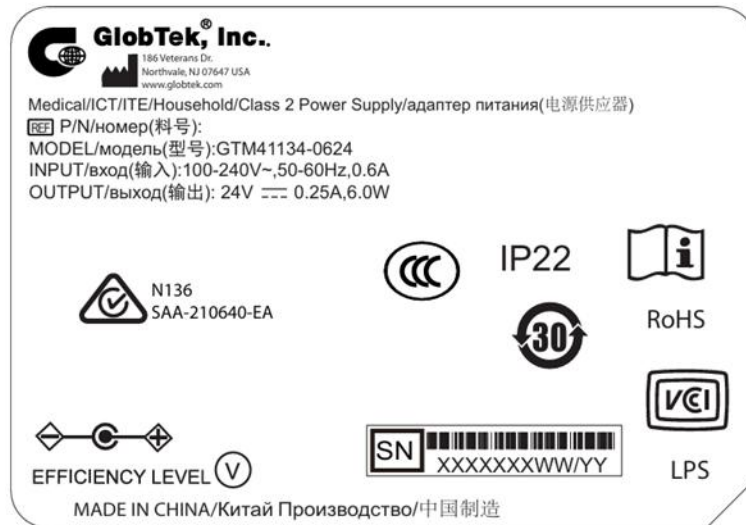

Rated output 7.5 W is only for model GTM96060-0706-1.0

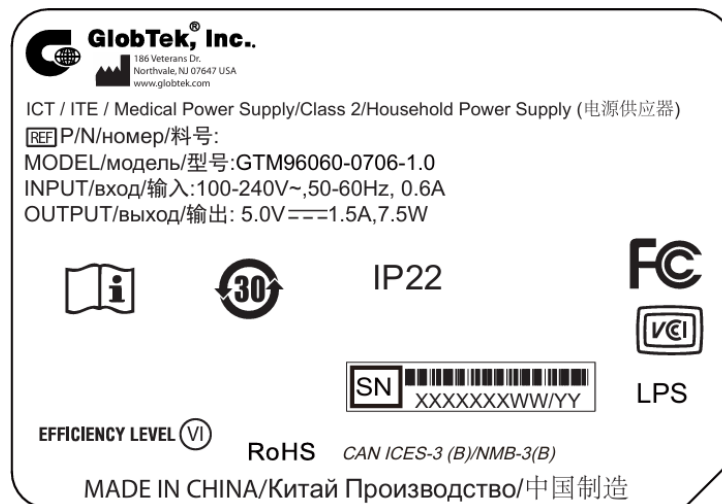

Note:

1, Markings of other models are similar as above except model name and ratings.



2, Symbols  and  for models with enclosure

3, 'YY' in SN represents manufacturing year; 'WW' in SN represents manufacturing week;

Manual/Service Instructions

- Operations and Service instructions are provided with the equipment. Doc name: _809-PSUQKITUSERMANUAL___ Revision: _Rev A6
See illustrations 5-12.

Alternate Listee Information

Alternate listees and product names or model numbers: None

Applicant's Responsibilities

Product Modifications:

Minor product modifications by the manufacturer may be allowed using the following guidelines:

1. Components identified in this report as "Listed, Recognized, or Certified" and **NOT** identified with a manufacturer name or part number may be exchanged with an alternate "Listed, Recognized, or Certified" component of equivalent value.

Example: Appliance Inlet Connector - Listed/Certified Component, IEC 320 style male connector, rated 250 volts and 20 amperes. Mechanically secured to the front panel with screws and locking washers.

- This inlet connector may be replaced with any Listed/Certified inlet connector with the same ratings as stated and where mechanical securement is maintained.

2. Components identified by a manufacturer name, part number, or with specific comments, (such as AC only, indoor use only, approved for use in this product only), may **NOT** be replaced or modified without prior approval from MET Laboratories.

Example: Circuit Breaker - Recognized/Certified Component, ABCD Co. P/N XYZ123, rated 250 volts maximum, 50/60 Hz, 25 full-load amperes, 31.3 trip amperes. Toggle handle marked with IEC on/off symbols. Mechanically secured to the front panel with screws and locking washers.

- This circuit breaker can **NOT** be modified or changed without prior approval by Eurofins E&E NA, LLC

This product was tested for use in Hazardous Locations. The manufacturer may not make any modifications without prior approval of Eurofins E&E NA, LLC

Applicant's Responsibilities (Continued)

Project Amendments:

For any changes related to product construction, manufacturing locations, if the product is intended to be marketed/sold under an alternate name or model number than that originally listed, or any issues which would require notification or change in the status of this file, please complete the form and return to Eurofins E&E NA, LLC following the instructions provided on the form.

For your convenience a Project Amendment Request (PAR) form is available for download at <http://corp.metlabs.com/safetyreq/> Alternatively, please provide it to your local Eurofins office or Eurofins Partner Representative.

If you are terminating or temporarily suspending production of this product for an extended period, please send a letter on company letterhead to:

Eurofins E&E NA, LLC
Attn: Follow Up Services Department
914 West Patapsco Avenue
Baltimore, Maryland 21230
USA
Fax: (410) 354-3313

Applicant's Responsibilities (Continued)

Manufacturing and Production-Line Tests and Documentation performed by Manufacturer.

All certified products are required to be subjected to production line testing as indicated below:

Dielectric Voltage-Withstand Test:

Each complete end product shall be capable of withstanding, without electrical breakdown, the application of a continuous sinusoidal or direct current voltage between uninsulated live parts and accessible dead metal parts that are likely to become energized in accordance with the following method.

Circuit Rating	Component Tested	Circuit Tested			
			Voltage (VAC)	Voltage (VDC)	Time (sec)
100-240 V	Main unit	Mains (L/N) to plug surface	2700	--	1
100-240 V	Main unit	Mains (L/N) to equipment surface	3600	--	1
100-240 V	Main unit	Mains (L/N) to output	3600	--	1

Documentation:

The manufacturer is required to record the production line test results. The data recorded is to include the type of test, date of test, serial number of the product, indications of pass, fail, or retest, test equipment utilized, calibration date of test equipment utilized, and the initials or signature of the test technician. Test records shall be required to be maintained from factory follow-up audit to factory follow-up audit and must be available for the inspectors' review. Records may be in the form of travelers, logs, computer files, or other such suitable documentation method.

Critical Components

Figure /item No.	Object/ Parts No.	Manufacturer / Trademark	Type/ Model	Technical Data	Standard (Edition / year)	Mark(s) of Conformity	Secured Method
1,2 /1	Plastic enclosure of plug and main unit	SABIC JAPAN LLC	945(GG)	PC, V-0, Min. thickness: 2.0mm, 120°C	UL 94 UL 746 A/B	UL/CSA Recognized (E207780)	Complete enclosure is fixed by ultrasonic welding. The plug was secured to main unit by snap fits.
	Alternate	SABIC INNOVATIVE PLASTICS B V	945(GG)	PC, V-0, Min. thickness: 2.0mm, 120°C	UL 94 UL 746 A/B	UL/CSA Recognized (E45329)	See Illustration 1 for details
-	Description	Enclosure: GTM96060:74 (L)*43.5(W)*36.8(H)+/- 1.0 mm GTM41134: 74 (L)*43.5(W)*35.3(H)+/- 1.0 mm					
3-7/2	PCB	WALEX ELECTRONIC (WUXI) CO LTD	T2, T2A, T2B, T4	min. V-0, 130°C	UL 796	UL/CSA Recognized (E154355)	Covered by plastic enclosure
	Alternate	Various	Various	min. V-0, 130°C	UL 796	UL/CSA Recognized	
3-7/3	Fuse (F1, F2) (F2 is optional)	SUZHOU WALTER ELECTRONIC CO LTD	2010	T1A, 250V	UL 248- 1 UL 248- 14	UL/CSA Recognized (E56092)	Soldered on PCB
	Alternate	Conquer Electronics Co., Ltd.	MST	T1A, 250V	UL 248- 1 UL 248- 14	UL/CSA Recognized (E82836)	
	Alternate	Dongguan Better	932	T1A, 250V	UL 248- 1 UL 248- 14	UL/CSA Recognized (E300003)	
	Alternate	Bel Fuse Ltd.	RST- Serie(s)	T1A, 250V	UL 248- 1 UL 248- 14	UL/CSA Recognized (E20624)	

Critical Components (Continued)

Figure /item No.	Object/ Parts No.	Manufacturer/ Trademark	Type/ Model	Technical Data	Standard (Edition / year)	Mark(s) of Conformity	Secured Method
7/4	Appliance inlet Class II units for FWT2 model	Zhejiang LECI Electronics Co., Ltd.	DB-8	2.5A, 250Vac, C8	UL 498	UL/CSA Recognized (E302229)	Intended to be fixed in end-product
	Alternate	Rich Bay Co., Ltd.	R-201SN90	2.5A, 250Vac, C8	UL 498	UL/CSA Recognized (E184638)	
	Alternate	Sun Fair Electric Wire & Cable (HK) Co. Ltd.	S-01	2.5A, 250Vac, C8	UL 498	UL/CSA Recognized (E226643)	
	Alternate	Rong Feng Industrial Co., Ltd.	RF-180	2.5A, 250Vac, C8	UL 498	UL/CSA L Recognized (E102641)	
	Alternate	Zhe Jiang Bei Er Jia Electronic Co., Ltd.	ST-A03 series	2.5A, 250Vac, C8	UL 498	UL/CSA Recognized (E225980)	
7/5	Insulating tube used on appliance inlet	SHENZHEN WOER HEAT-SHRINKABLE MATERIAL CO LTD	RSFR RSFR-H RSFR-HPF	600V, 125°C, VW-1	UL 224	UL/CSA Recognized (E203950)	Fixed on internal wire
	Alternate	QIFURUI ELECTRONICS CO	QFR-h	600V, 125°C, VW-1	UL 224	UL/CSA Recognized (E225897)	
	Alternate	DONGGUAN SALIPT CO LTD	SALIPT S-901-300 SALIPT S-901-600	Min. 300V, 125°C, VW-1	UL 224	UL/CSA Recognized (E209436)	
	Alternate	GUANGZHOU KAIHENG ENTERPRISE GROUP	K-2 (+) K-2 (CB)	Min. 300V, 125°C, VW-1	UL 224	UL/CSA Recognized (E214175)	
	Alternate	CHANGYUAN ELECTRONICS (SHENZHEN) CO LTD	CB-HFT	Min. 300V, 125°C, VW-1	UL 224	UL/CSA Recognized (E180908)	
	Alternate	SHENZHEN WOLIDA TRADING CO LTD	RSFR-H	600V, 125°C, VW-1	UL 224	UL/CSA Recognized (E329530)	

Critical Components (Continued)

Figure/ item No.	Object/ Parts No.	Manufacturer/ Trademark	Type/ Model	Technical Data	Standard (Edition / year)	Mark(s) of Conformity	Secured Method
3-7/6	Bridging capacitor (CY1, CY2) (Optional)	TDK Corporation	CD	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E37861)	Soldered on PCB
	Alternate	Success Electronics Co., Ltd.	SE	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E114280)	
	Alternate	Success Electronics Co., Ltd.	SB	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E114280)	
	Alternate	Murata Mfg. Co., Ltd.	KX	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E37921)	
	Alternate	SHANTOU HIGH-NEW TECHNOLOG Y DEVELOPMN T ZONE SONGTIAN ENTERPRISE CO LTD	CD	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E208107)	
	Alternate	JYA-NAY Co., Ltd.	JN	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E201384)	
	Alternate	JUSUN (TAISHAN) ELECTRONIC S LTD	JB-series	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E253194)	
	Alternate	XIANGTAI ELECTRONIC (SHENZHEN) CO LTD	YO series	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E319473)	
	Alternate	Dongguan Easy-gather Electronic Co., Ltd.	DCF	Y1, min. 250VAC, 125°C max. 4700pF	UL 1414	UL/CSA Recognized (E252221)	

Critical Components (Continued)

Figure/ item No.	Object/ Parts No.	Manufacturer/ Trademark	Type/ Model	Technical Data	Standard (Edition / year)	Mark(s) of Conformity	Secured Method
3-7/7	Varistor (MOV1) (optional)	Xiamen Set Electronics Co., Ltd.	TFV8S471 K TFV10S47 1K	Max. Continuous voltage: Min. 300Vac(rms), Min. 105°C, The coating is Min. V-0	UL 1449	UL/CSA Recognized (E322662)	Soldered on PCB
	Alternate	SHANTOU HIGH-NEW TECHNOLOG Y DEVELOPMN T ZONE SONGTIAN ENTERPRISE CO LTD	10D621K 14D621K	Max. Continuous voltage: Min. 385Vac(rms), Min. 125°C, The coating is Min. V-0	UL 1449	UL/CSA Recognized (E330837)	
	Alternate	Guangdong Huiwan Electronics Technology Co.Ltd.	V-621K-10 DEH V-621K-14 DEH	Max. Continuous voltage: Min. 385Vac(rms), Min. 125°C, The coating is Min. V-0	UL 1449	UL/CSA Recognized (E480104)	
	Alternate	Thinking Electronic Industrial Co., Ltd.	TVR10621 TVR14621	Max. Continuous voltage: Min. 385Vac(rms), Min. 125°C, The coating is Min. V-0	UL 1449	UL/CSA Recognized (E314979)	

Critical Components (Continued)

Figure/ item No.	Object/ Parts No.	Manufacturer/ Trademark	Type/ Model	Technical Data	Standard (Edition / year)	Mark(s) of Conformity	Secured Method
1,2/8	Output cord	Various	Various	Min. 22AWG, min. 300Vac, min. 80°C Min. VW-1	UL 758	UL/CSA	Fixed by enclosure
8-10/9	Transfor mer (T1)	GlobTek/ENG BOAM / HAOPUWEI	--	Class B	AAMI ES 60601-1, CSA C22.2 No. 60601- 1	Tested with equipment	Soldered on PCB
8-10 /10a	- Insulation system	GLOBTEK INC	GTX-130- TM	Class B	UL 1446	UL/CSA Recognized (E243347)	Fixed in transformer
	Alternate	ENG Electric Co Ltd	ENG130-1	Class B	UL 1446	UL/CSA Recognized (E308897)	
	Alternate	HAOPUWEI	ZT-130	Class B	UL 1446	UL/CSA Recognized (E315575)	
	Alternate	BOAM	BOAM-01, or B1	Class B	UL 1446	UL/CSA Recognized (E252379)	
8-10 /10b	-Magnet wire (primary)	WUXI JUFENG COMPOUND LINE CO LTD	2UEWB	MW75#, 130°C	UL 1446	UL/CSA Recognized (E206882)	Fixed in transformer
	Alternate	JIANGSU DARTONG M & E CO LTD	UEW	MW 75-C, 130°C	UL 1446	UL/CSA Recognized (E237377)	
	Alternate	SHANDONG SAINT ELECTRIC CO LTD	UEW/130	MW75#, 130°C	UL 1446	UL/CSA Recognized (E194410)	
	Alternate	NINGBO JINTIAN NEW MATERIAL CO LTD	2UEW	MW 79#, 155°C	UL 1446	UL/CSA Recognized (E227047)	
	Alternate	Various	Various	Min 130 °C	UL 1446	UL/CSA Recognized	

Critical Components (Continued)

Figure/ item No.	Object/ Parts No.	Manufacturer/ Trademark	Type/ Model	Technical Data	Standard (Edition / year)	Mark(s) of Conformity	Secured Method
8-10 /10c	-Triple Insulation wire	GREAT LEOFLON INDUSTRIAL CO LTD	TRW(B) series	130°C	UL 2353	UL/CSA Recognized (E211989)	Fixed in transformer
	Alternate	TOTOKU INC.	TIW-2X\$+	130°C	UL 2353	UL/CSA Recognized (E166483)	
	Alternate	HOI LUEN ELECTRICAL MFR CO LTD	THL-F-xx, THL-F-SB- xx	130°C	UL 2353	UL/CSA Recognized (E257525)	
	Alternate	Furukawa Electric Co., Ltd.	TEX-E	130°C	UL 2353	UL/CSA Recognized (E006735)	
8-10 /10d	-Bobbin	CHANG CHUN PLASTICS CO LTD	T375J, T375HF	V-0, 150°C, thickness 0.71 mm min.	UL 796	UL/CSA Recognized (E59481)	Fixed in transformer
	Alternate	Resonac Techno Service Corporation	CP-J-8800	Phenolic, V-0, Min. 150°C, Min. 0.71mm thickness	UL 796	UL/CSA Recognized (E514814)	
	Alternate	SUMITOMO BAKELITE CO LTD	PM-9820, PM-9630	V-0, 150°C, thickness 071 mm min.	UL 796	UL/CSA Recognized (E41429)	

Critical Components (Continued)

Figure/ item No.	Object/ Parts No.	Manufacturer/ Trademark	Type/ Model	Technical Data	Standard (Edition / year)	Mark(s) of Conformity	Secured Method
8-10 /10e	-Tape	3M COMPANY	1350F-1, 1350T-1 44	130°C	UL 510	UL/CSA Recognized (E17385)	Fixed in transformer
	Alt.	BONDTEC PACIFIC CO LTD	370S	130°C	UL 510	UL/CSA Recognized (E175868)	
	Alt.	JINGJIANG YAHUA PRESSURE SENSITIVE GLUE CO LTD	PZ CT WF	130°C	UL 510	UL/CSA Recognized (E165111)	
	Alt.	JINGJIANG JINGYI ADHESIVE PRODUCT CO LTD	JY25-A	130°C	UL 510	UL/CSA Recognized (E246950)	
	Alt.	CHANG SHU LIANG YI TAPE INDUSTRY CO LTD	LY-XX	130°C	UL 510	UL/CSA Recognized (E246820)	
8-10 /10f	-PTFE tubing	Great Holding Industrial Co Ltd	TFT, TFS	Min. 300V, 200°C	UL 224	UL/CSA Recognized (E156256)	Fixed in transformer
	Alt.	Shenzhen Woer Heat- Shrinkable Material Co Ltd	WF	600V, 200°C	UL 224	UL/CSA Recognized (E203950)	
	Alt.	Changyuan Electronics (Shenzhen) Co Ltd	CB-TT-T, CB-TT-S	Min. 300V, 200°C	UL 224	UL/CSA Recognized (E180908)	
1,2/11	Integral plug	GlobTek, Inc.	--	125V, 1-15P	UL 498 NEMA 1- 15	Intertek (220201762 SHA-001)	Clapped in plastic enclosure of main unit


Critical Drawings

Title:	Drawing No.:	Rev. Level:	Date:
Circuit diagram GTM96060	Illustration 1	B	2019-02-12
PCB layout GTM96060 EL6	Illustration 2	A8	2025-03-19
Circuit diagram GTM41134	Illustration 3	03	2020-09-02
PCB layout GTM41134	Illustration 4	06	2018-12-04
Transformer	Illustration 13-14	A.0	2017-04-21

Figures

Figure 1.


External view of GTM96060-0706-1.0, models with enclosure in series GT*41134***** and GT*96060*****

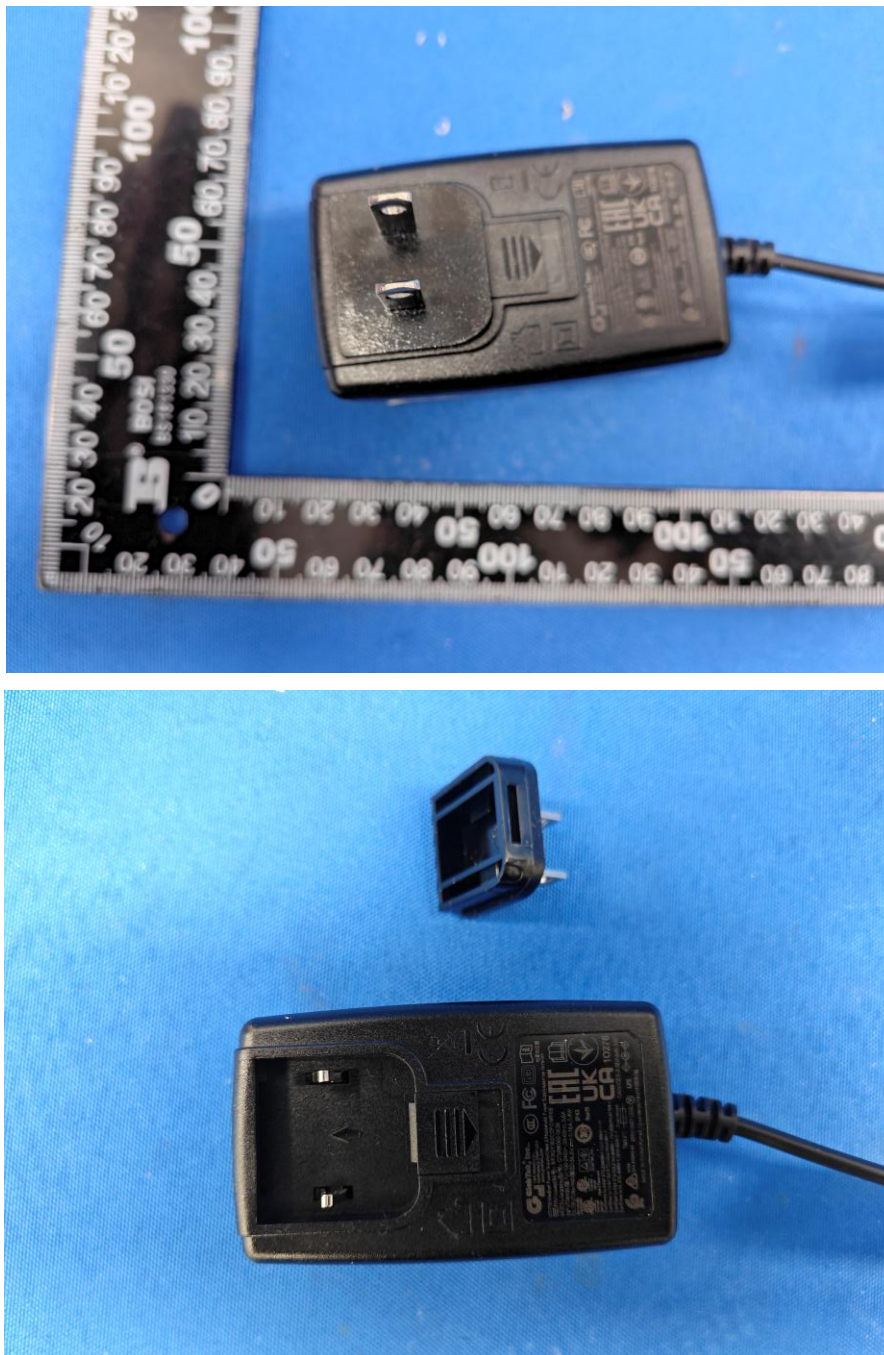

Figures (Continued)

Figure 2.

External view of GTM96060-0706-1.0, models with enclosure in series GT*41134***** and GT*96060*****

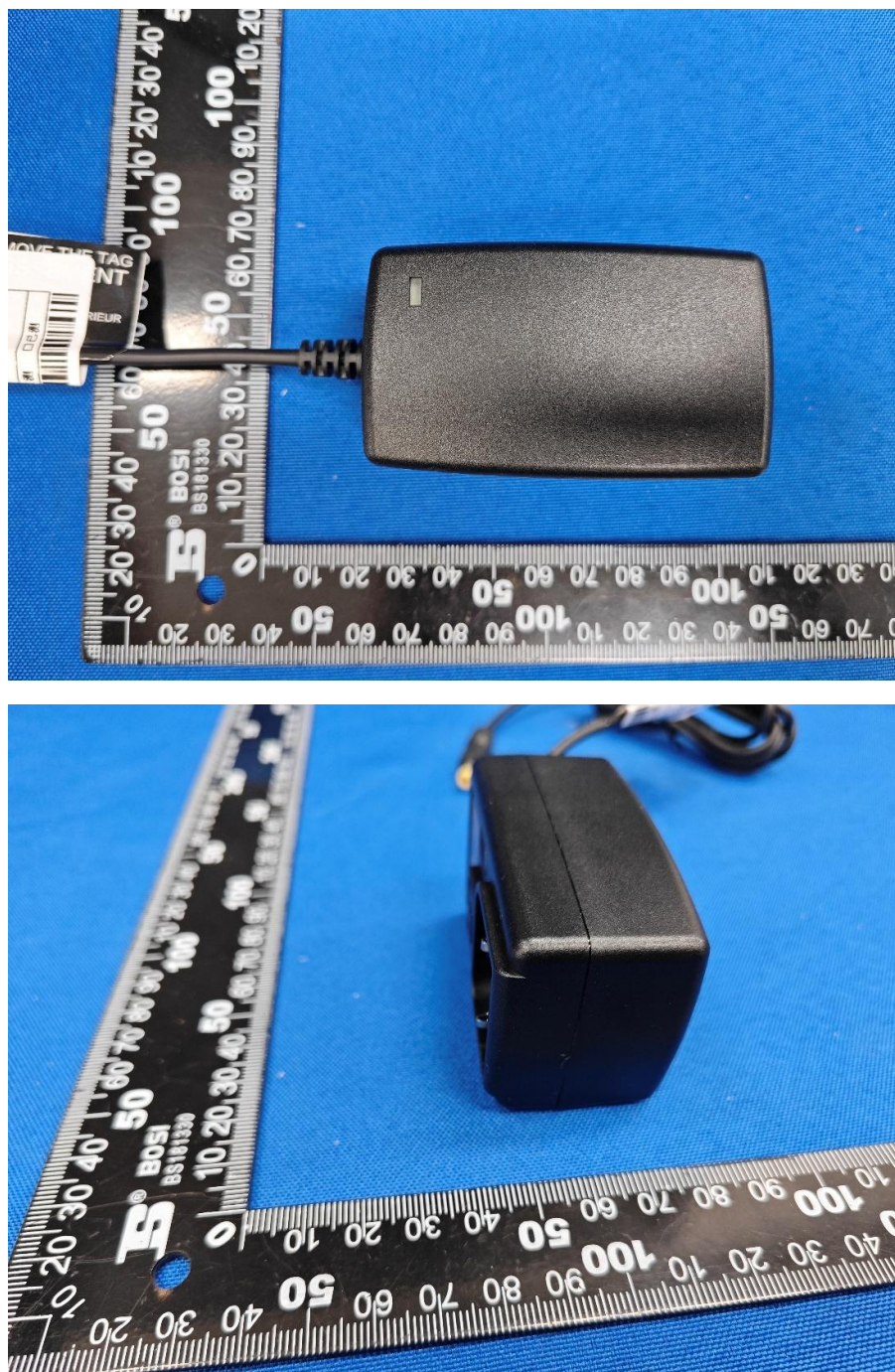


Figures (Continued)

Figure 3.

Internal view of GT*96060***** series, GTM96060-0706-1.0





Internal view of GT*41134***** series

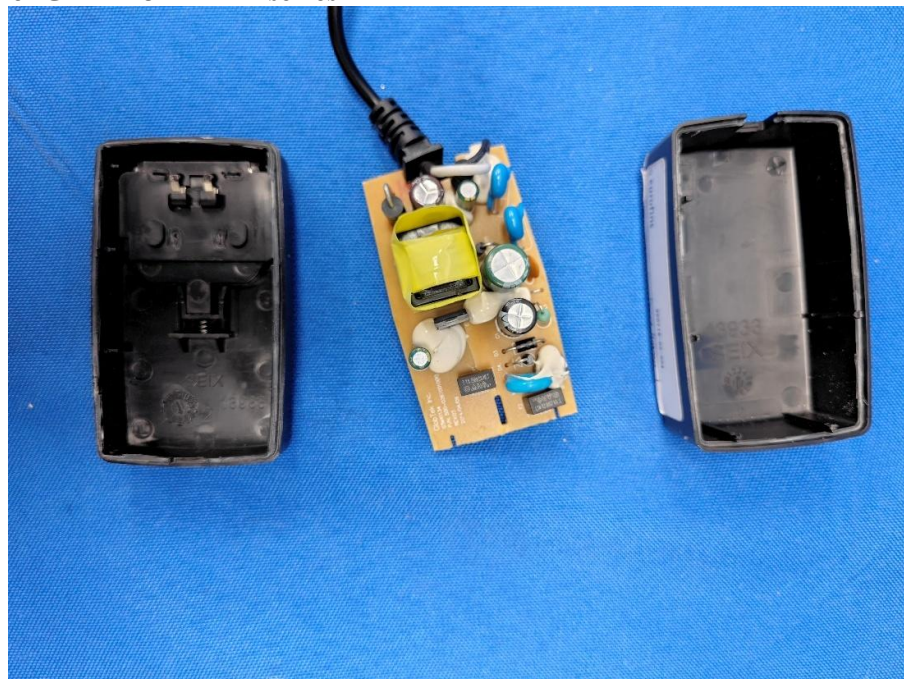

Figures (Continued)

Figure 4.



PCB of GT*96060***** series, GTM96060-0706-1.0

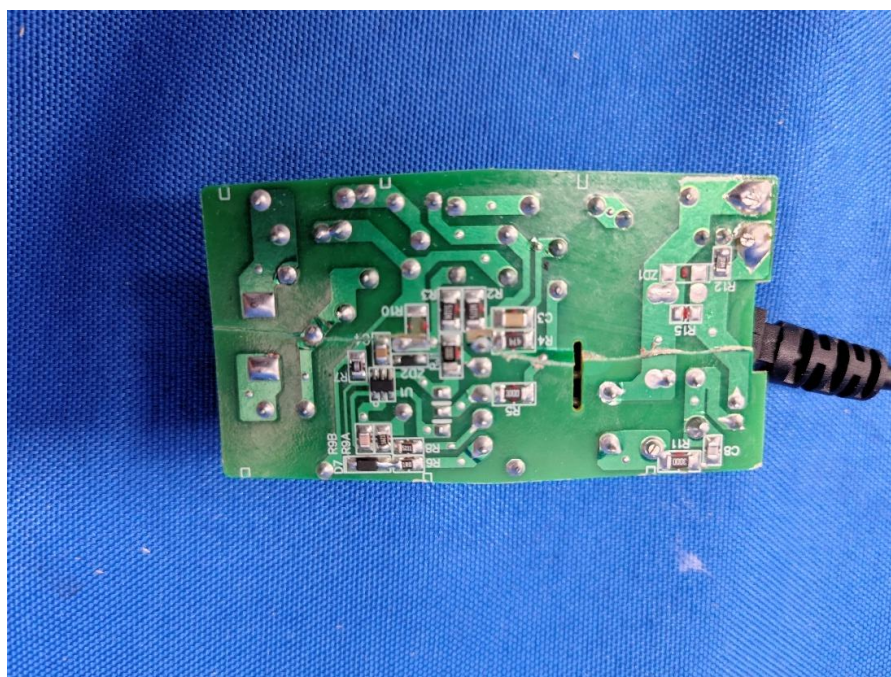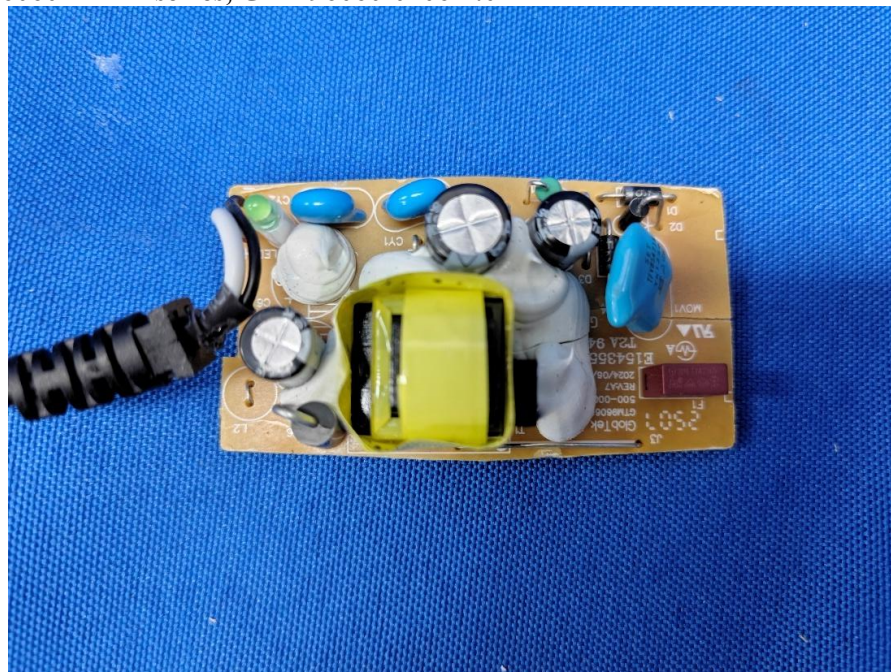

Figures (Continued)

Figure 5.


PCB of GT*41134 series

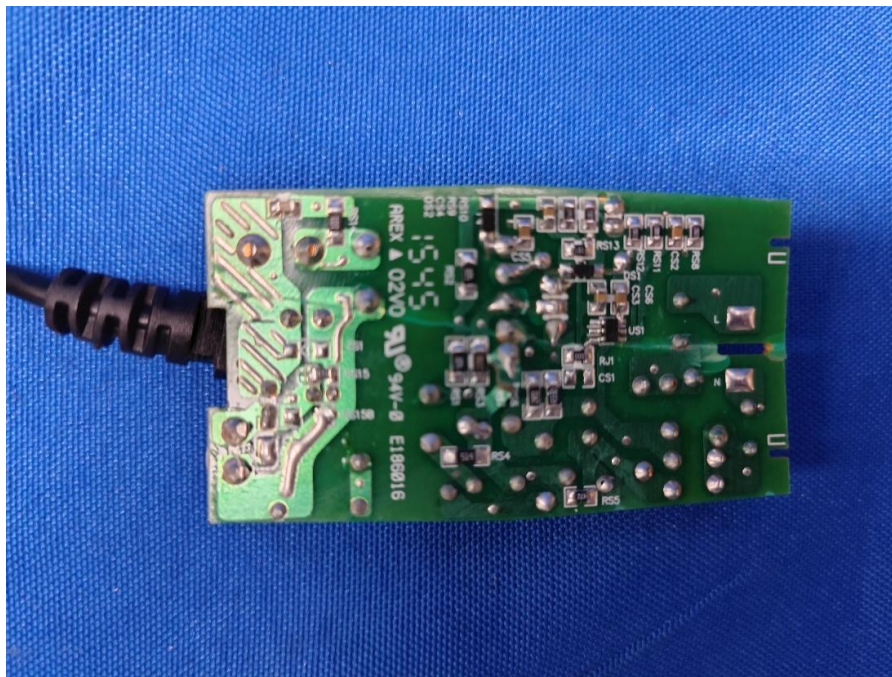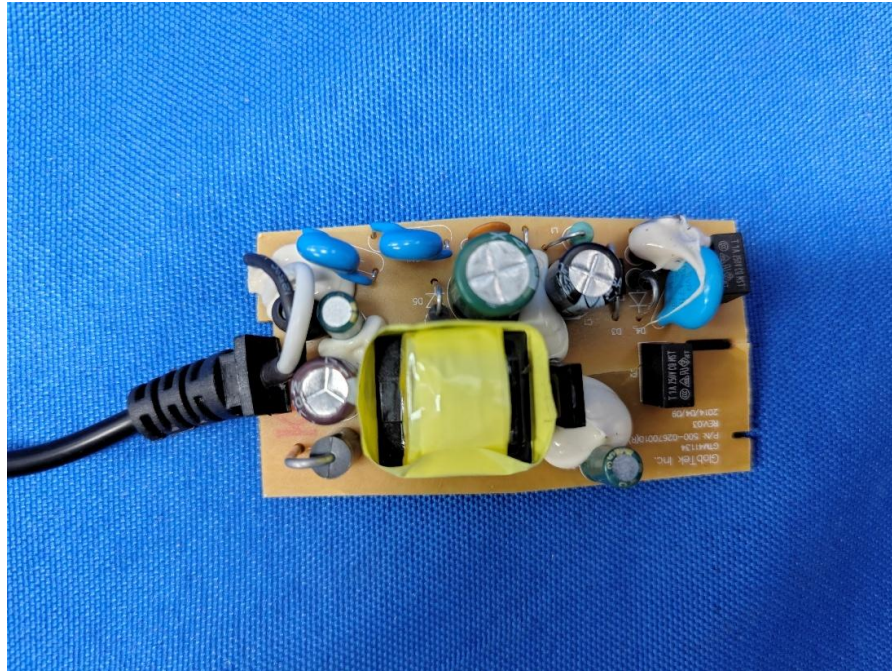

Figures (Continued)

Figure 6.



Open frame models with connector: GT*41134****-FW*, GT*96060****-FW*

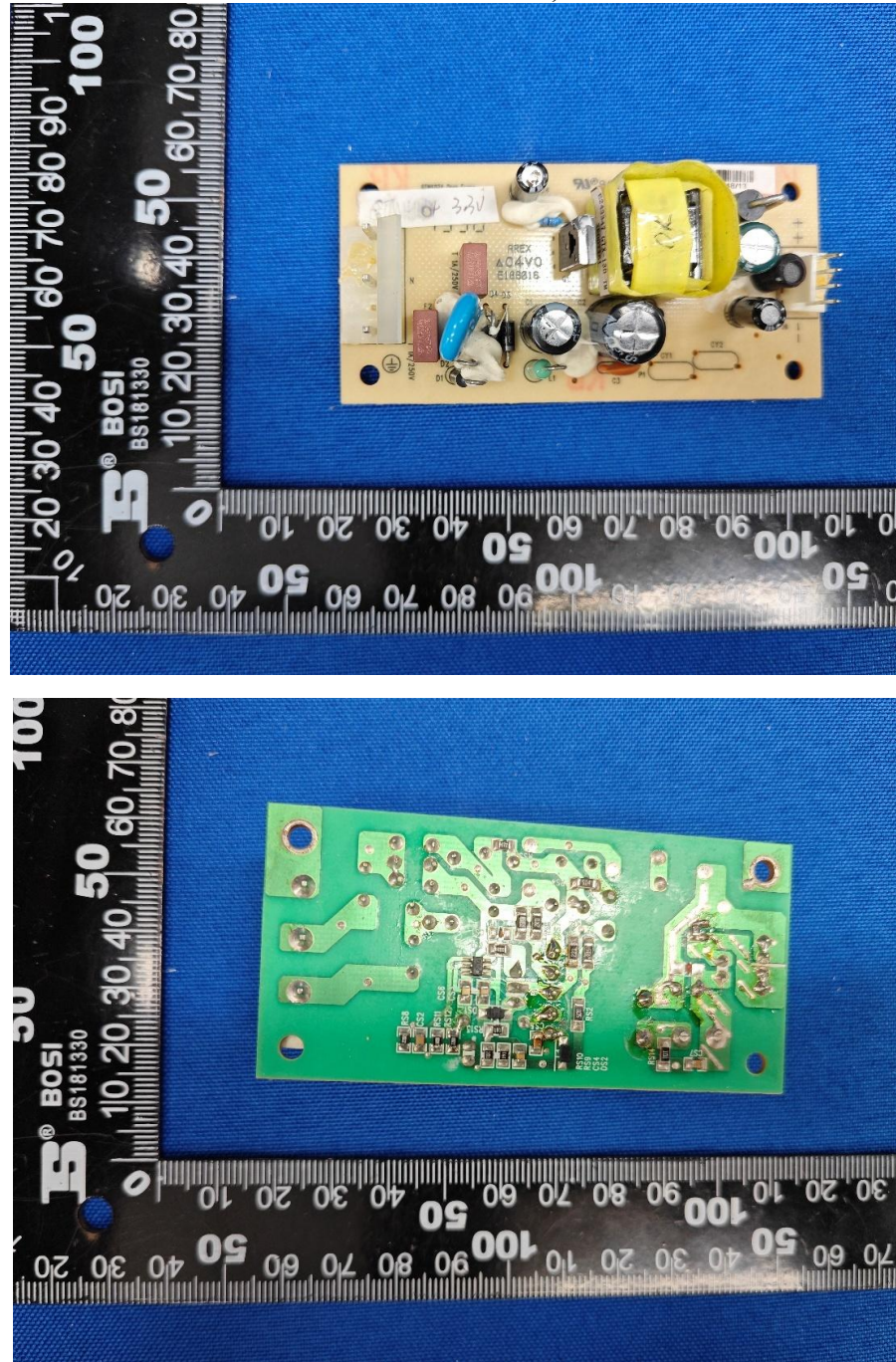

Figures (Continued)

Figure 7.



Open frame models with appliance inlet GT*41134****-FWT2*, GT*96060****-FWT2*

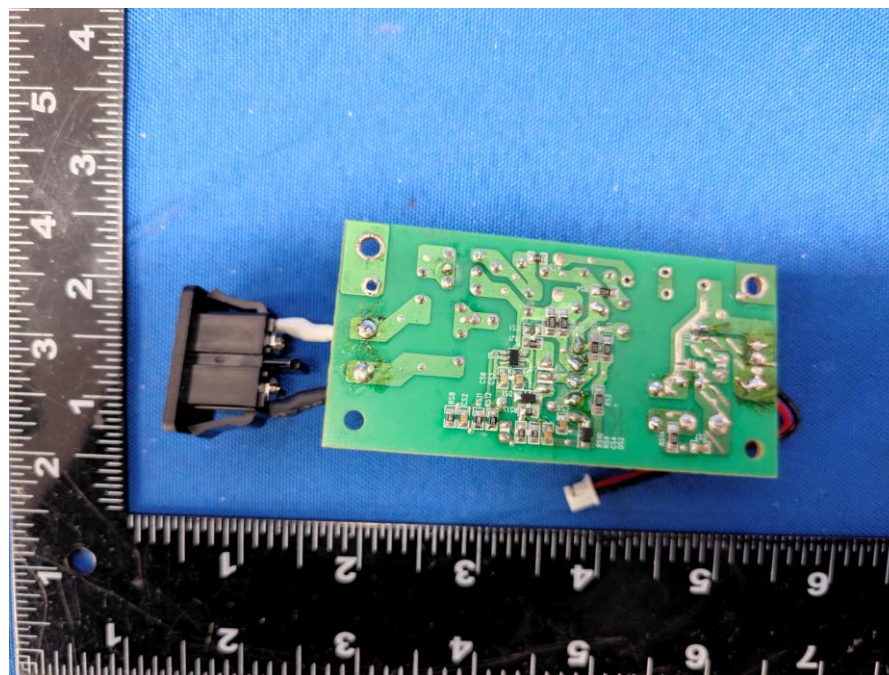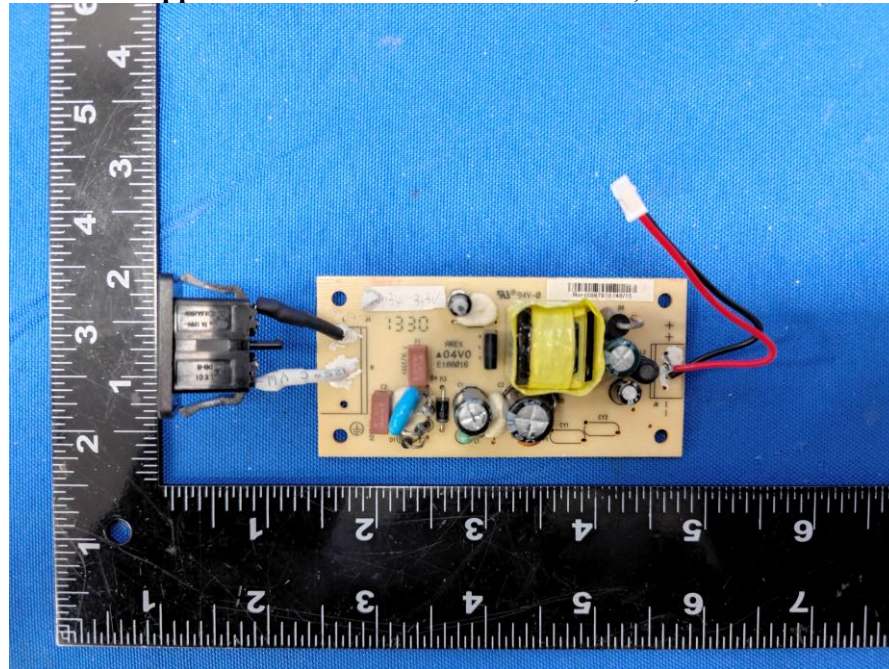

Figures (Continued)

Figure 8.



Transformer

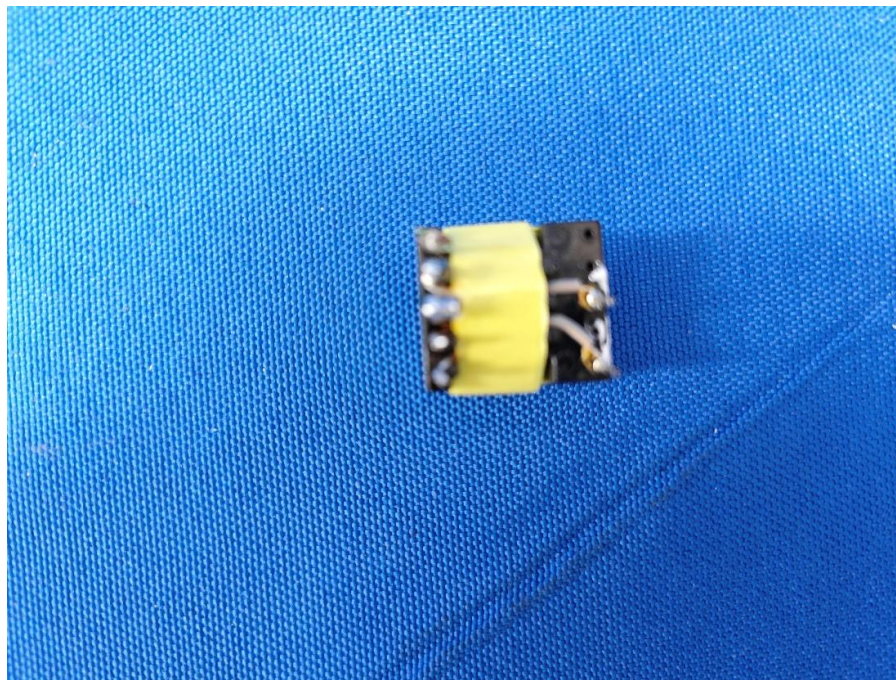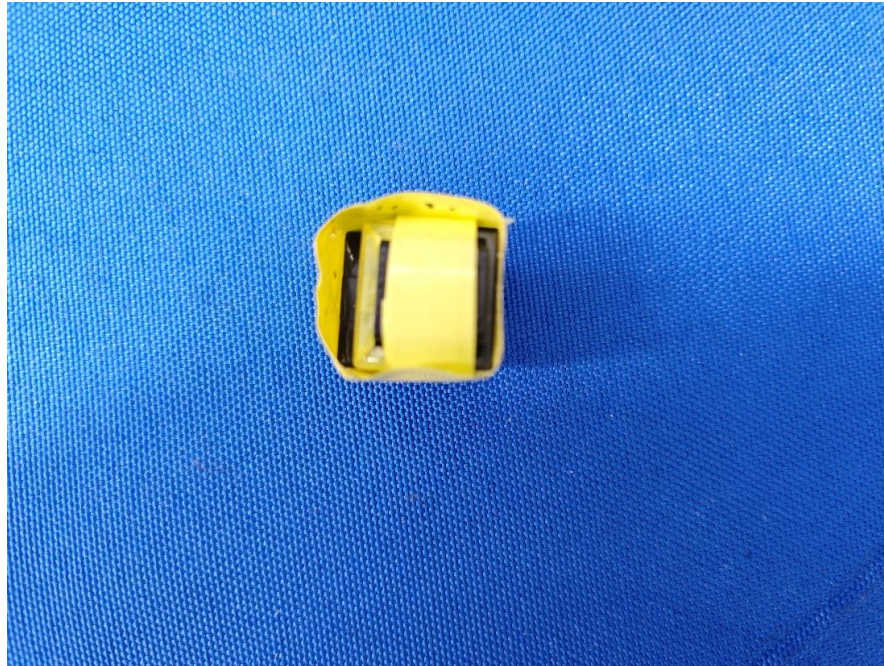

Figures (Continued)

Figure 9.

Transformer

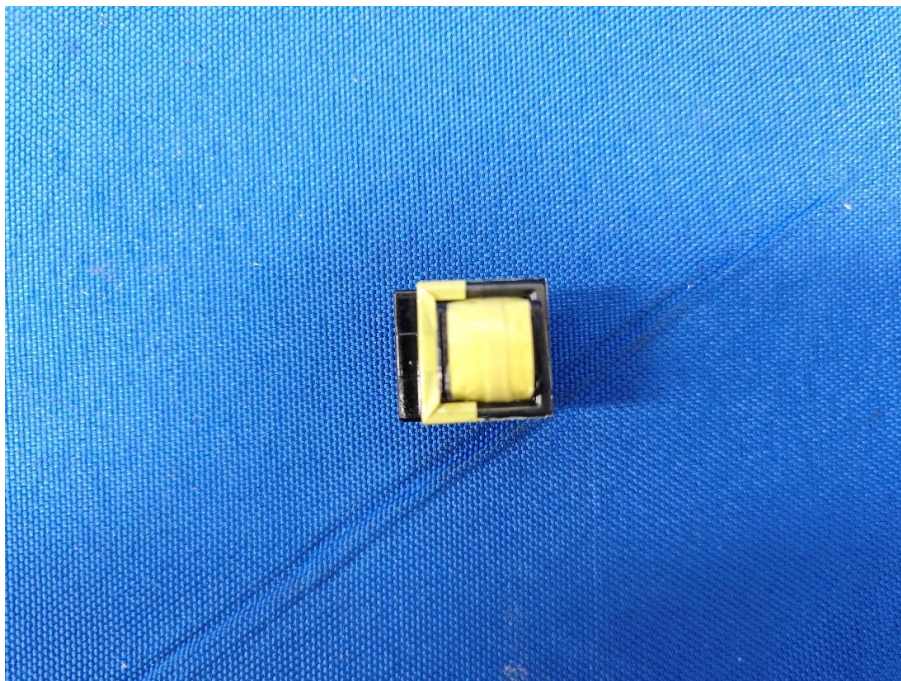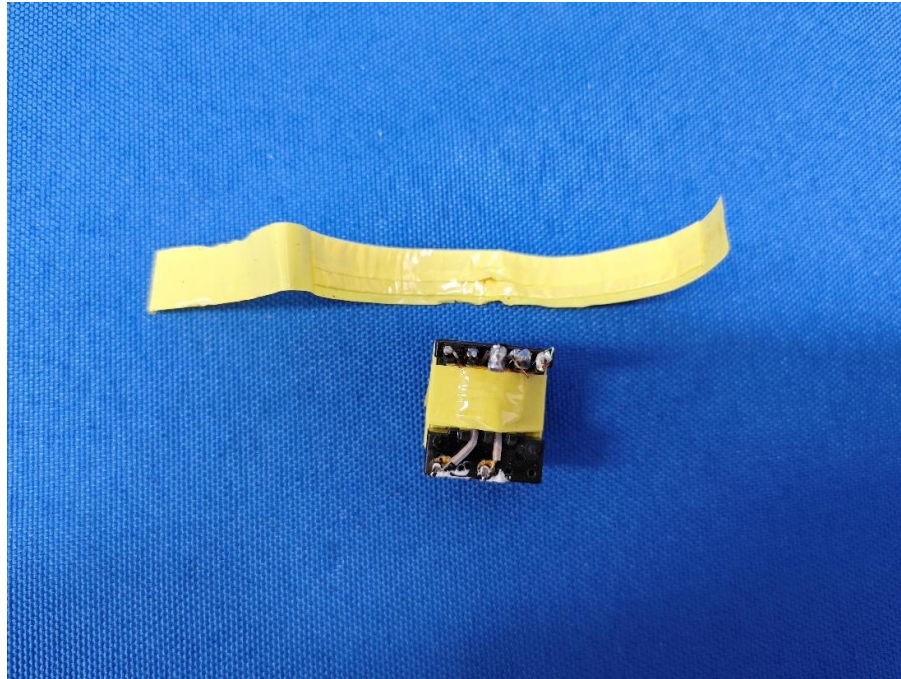

Figures (Continued)

Figure 8.


Transformer



Illustrations

Illustration 1.


Circuit diagram GTM96060

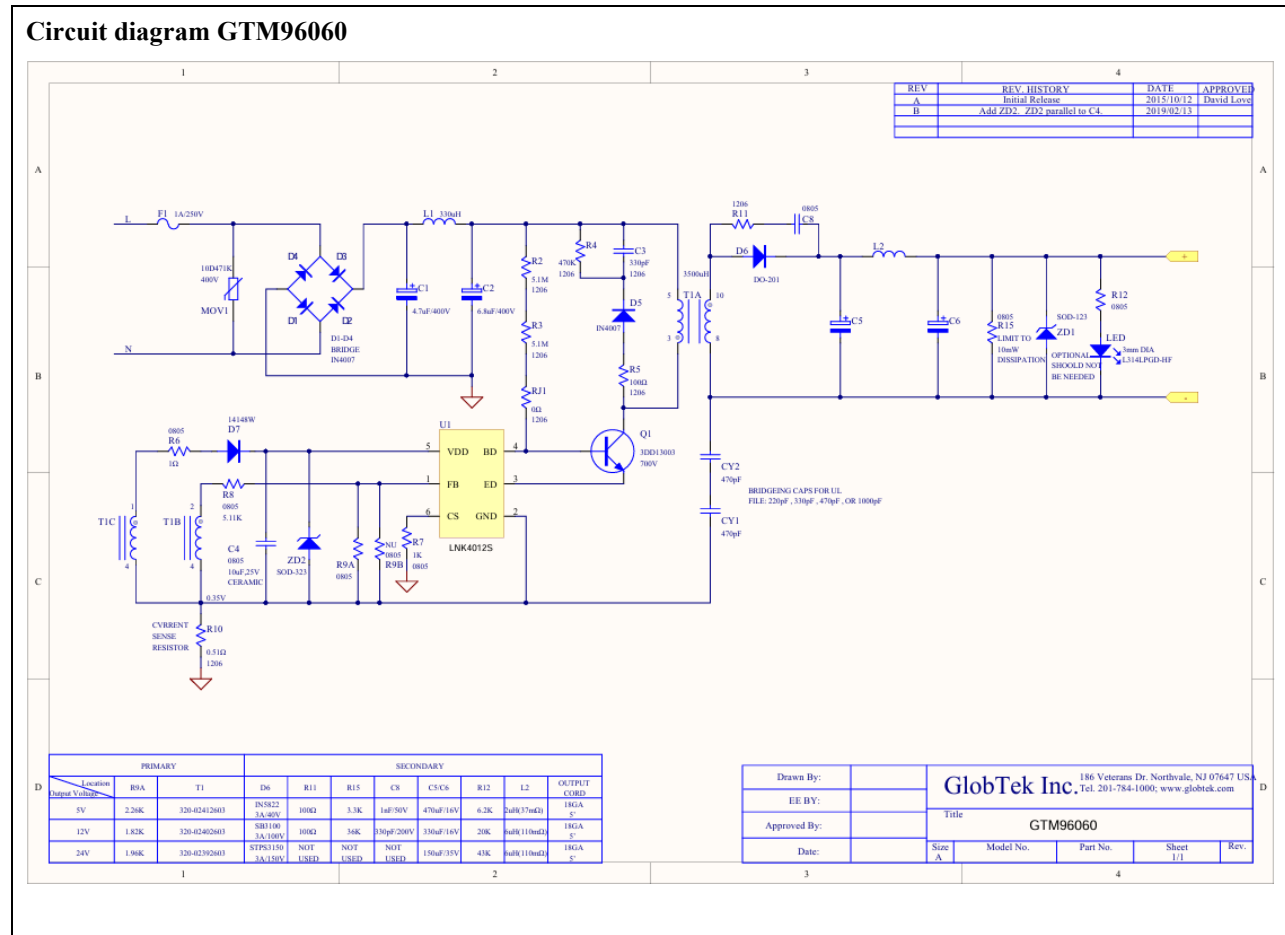


Illustrations (Continued)

Illustration 2.

PCB layout GTM96060

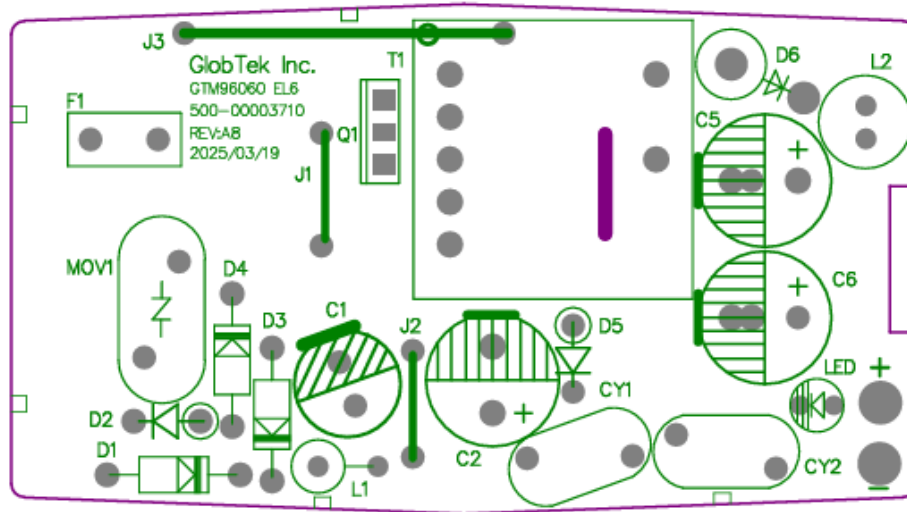


PN: 500-00003710(R) PCB MASK COLOR: Green
 Material: FR-1 1.6mm T Tolerance: +/- 0.1(mm)
 Copper: 1oz. REV: A8
 Top Text: White Flame Class: 94V-0
 Bottom Text: White Surface Dispose: OSP
 PCB VENDOR: ADD UL APPROVED MANUFACTURER LOGO,
 PCB MATERIAL & FLAMMABILITY & YOUR COMPANY LOGO.

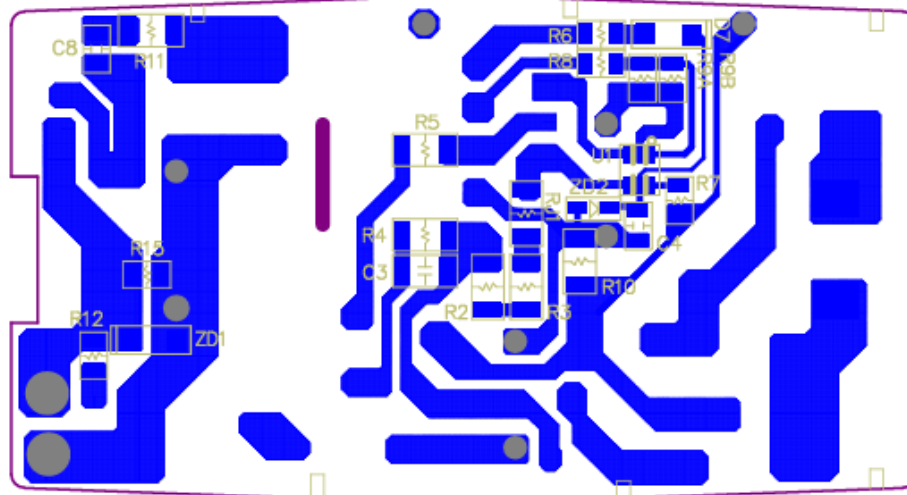


Illustrations (Continued)

Illustration 3.

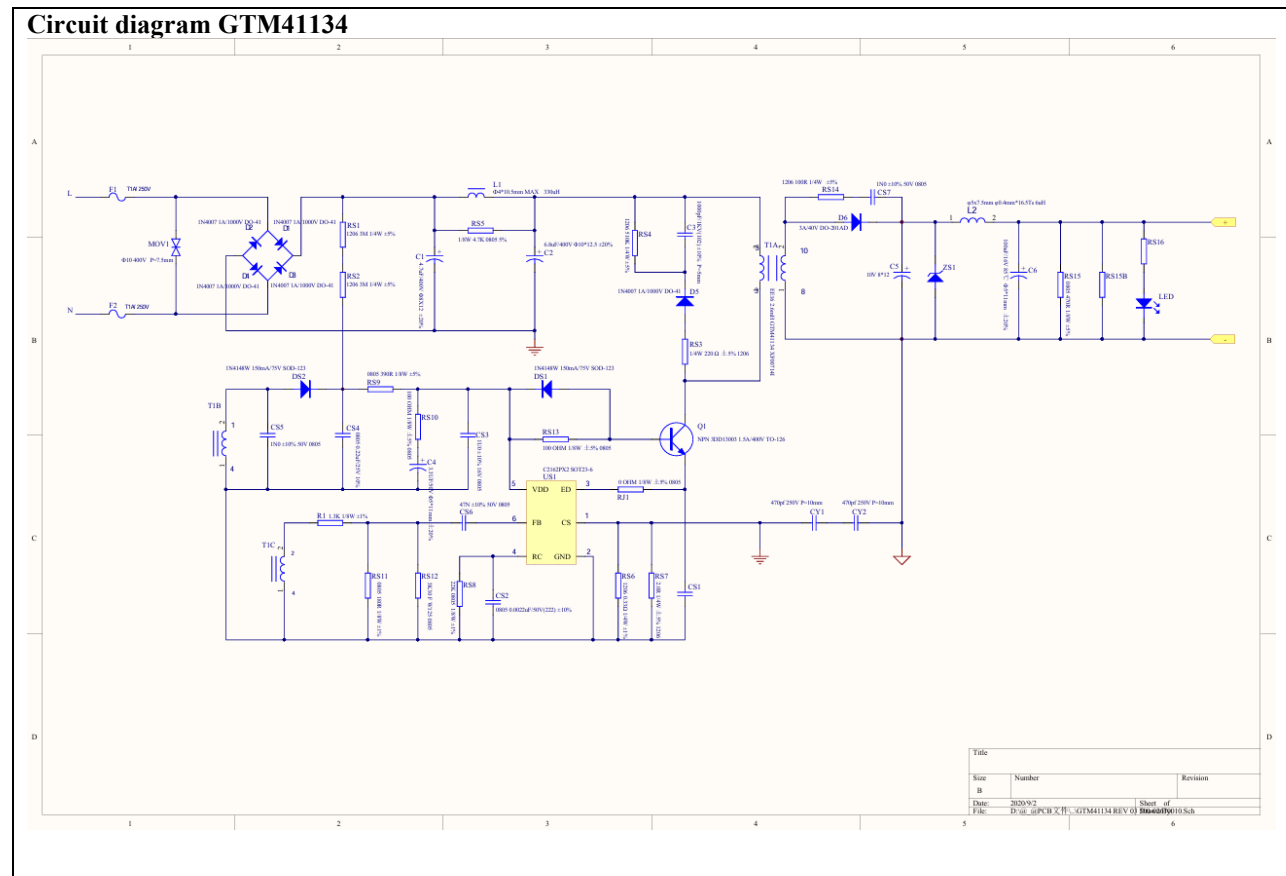

Illustrations (Continued)

Illustration 4.

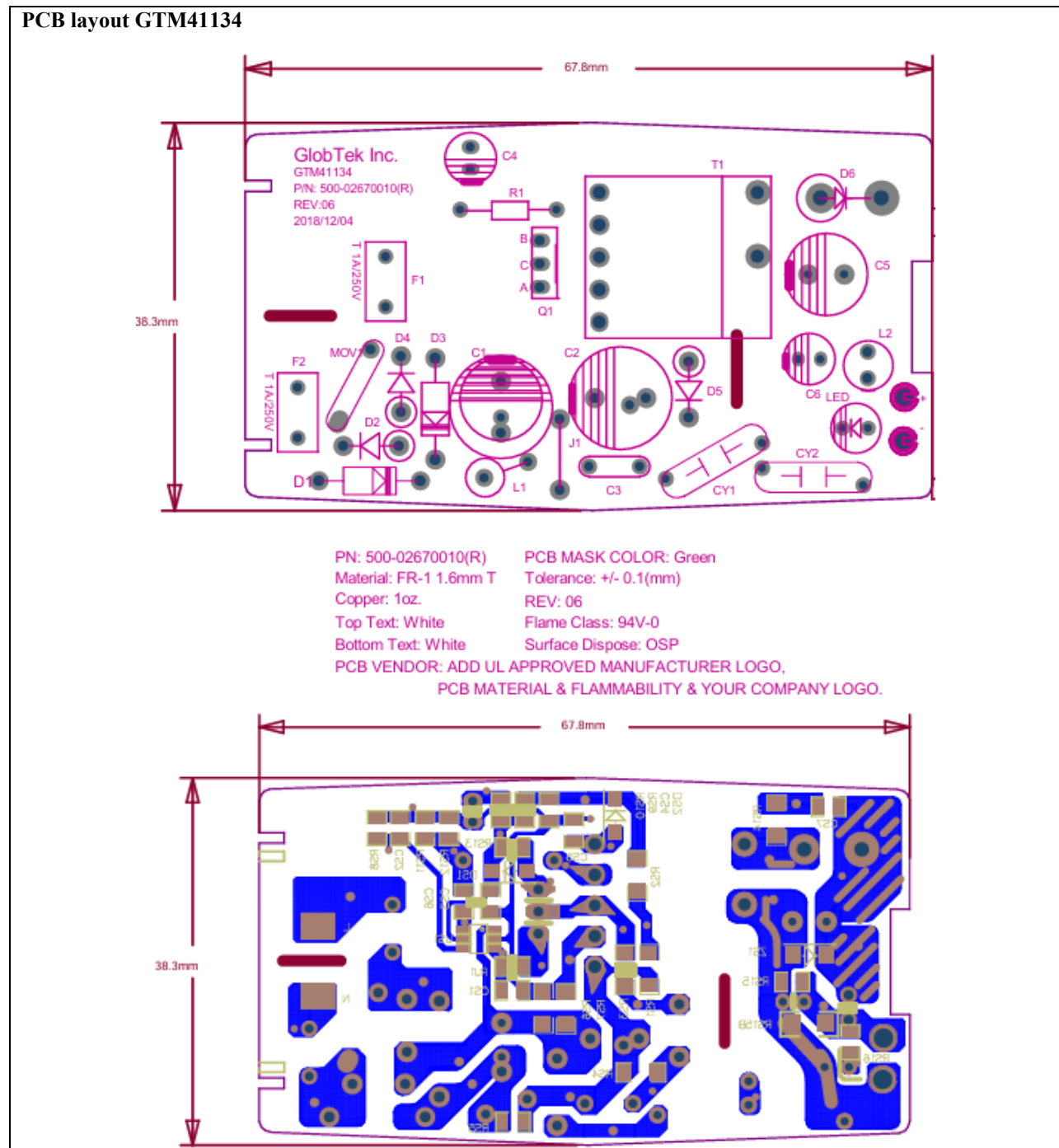

Illustrations (Continued)

Illustration 5.

User manual

GlobTek, Inc.

User manual for Model No.: GTM41134 & GTM96060 series

ELECTRICAL SPECIFICATIONS:

Input: 100-240Vac 50-60Hz 0.3A or 0.6A

Output: 3.3V-36.0Vdc; Max. 1.8A; Max. 7.5W, details refer to making label.

MTBF: 1,500,000 HOURS AT 25°C ambient temperature

Operating Temperature: -10° C TO 50 ° C

Storage Temperature: -25° C TO 70 ° C

Humidity: 0% TO 95% relative humidity

Ambient pressure range: 62-106kPaWeight: 0.13kg

Enclosure Size: 74 x 43.5 x 36.8 mm +/-1.0mm

Input connector: AC power receptacle with 2-Pin Plug

Output cord: No. 24 AWG min., VW-1, min.80 degree C, min.30 V;

1. Please read these safety instructions carefully.
2. This unit is for indoor use only, is expected to be used in IT equipment or Audio and Video products.
3. Please refer to the marking label on the unit for input and output ratings. Do not overload the power supply.
4. The socket-out shall be installed near the equipment and shall be easily accessible.
5. A readily accessible disconnect device shall be incorporated in the fixed wiring.
6. Disconnect the appliance coupler to separate the unit from mains supply.
7. The apparatus shall not be exposed to dripping or splashing.
8. Expected service life of this equipment is 5 years.

9. Output polarity: 

10. Class II equipment: 

11. Indoor use: 

12. Attention, Consult Accompanying Documents as applicable) 

13. Do not throw this electronic device into the trash when discarding. To minimize pollution and ensure utmost protection of the global environment, please recycle it.



14. The power supply has to be disposed appropriately: Please refer to local regulations (Waste Electrical and Electronic Equipment).

15. The LED is operating indicator (LED is optional)

16. Regardless to the legal warranty the manufacturer grants in accordance with the laws of its country, at least 1 year (in Germany 2 years). Begin of warranty is the sales date of the equipment to the final consumer.

The warranty covers only defects/deficiencies, which are due to failures of material or manufacturing.

Illustrations (Continued)

Illustration 6.

GlobTek, Inc.

Warranty repairs may be accomplished exclusively by an authorized customer service. In order to make your warranty claim the original (sales voucher) receipt (with sales date) must be attached.

Excluded from warranty are:

- Usual abrasion
- Inappropriate applications, like for example overloading of the equipment, not certified accessories and modified accessories.
- Damages by foreign effects, use of force or by foreign objects
- Damage, which results from neglect of the user manual, for example connection to a wrong voltage or neglect of the assembly instructions
- Complete or partly dismantled devices

17. The apparatus is not intended to use in oxygen-enriched environment, anesthesia environment.

18. The apparatus shall be complied with IEC 60601-1-2 electromagnetic.

19. Do not wrap the cable around your neck to cause asphyxiation

20. If the equipment doesn't work, please check whether the connection is firm firstly, other further questions, sending to the warranty point is recommended

21. Unauthorized modification of the equipment may cause hazards.

Illustrations (Continued)

Illustration 7.

GlobTek, Inc.

Manuel d'utilisation pour le modèle non:
GTM41134 & GTM96060 Les séries

Spécifications électriques:
Entrée en ligne: 100-240Vac 50-60Hz 0.3A or 0.6A

La production: 3.3V-36.0Vdc; Max. 1.8A; Max. 7.5W, les détails se réfèrent à faire l'étiquette.
MTBF: 1,500,000 HOURS AT 25°C ambient temperature
Température de fonctionnement: -10° C TO 50 ° C
Température de stockage: -25° C TO 70 ° C
humidité: 0% TO 95% humidité relative
Plage de pression ambiante: 62-106kPa
Poids poids: 0.13kg
Taille du boîtier: 74 x 43.5 x 36.8 mm +/-1.0mm
Connecteur d'entrée: Prise de courant AC avec prise 2 broches
Cordon de sortie: No. 24 AWG min., VW-1, min.80 degree C, min.30 V;

1. Veuillez lire attentivement ces consignes de sécurité.
2. Cette unité est pour l'usage d'intérieur seulement, devrait être employée dans l'équipement informatique ou les produits Audio et vidéo.
3. Veuillez vous référer à l'étiquette de marquage sur l'unité pour les cotes d'entrée et de sortie. Ne surchargez pas l'alimentation.
4. La chaussure doit être installée près de l'équipement et doit être facilement accessible.
5. Un dispositif de déconnexion facilement accessible doit être incorporé dans le câblage fixe.
6. Débranchez le coupleur de l'appareil pour séparer l'unité de l'alimentation secteur.
7. L'appareil ne doit pas être exposé à l'égouttement ou aux éclaboussures.
8. La durée de vie prévue de cet équipement est de 5 ans.
9. Polarité de sortie: 
10. Équipement de classe II: 
11. Utilisation à l'intérieur: 
12. Attention, consulter les Documents d'accompagnement - le cas échéant
13. Ne jetez pas cet appareil électronique à la poubelle lors du rejet. Pour minimiser la pollution et assurer une protection maximale de l'environnement mondial, veuillez le recycler.
14. L'alimentation doit être éliminée de manière appropriée: veuillez vous référer à la réglementation locale (déchets d'équipements électriques et électroniques).
15. La LED fonctionne indicateur (la LED est facultative).
16. Indépendamment de la garantie légale que le fabricant accorde conformément aux lois de son pays, au moins 1 an (en





Illustrations (Continued)

Illustration 8.

GlobTek, Inc.

Allemagne 2 ans). Début de la garantie est la date de vente de l'équipement au consommateur final.

La garantie ne couvre que les défauts/déficiences, qui sont dus à des échecs de matériel ou de fabrication.

Les réparations sous garantie peuvent être effectuées exclusivement par un service à la clientèle autorisé. Afin de faire votre demande de garantie, le reçu original (bon de vente) (avec la date de vente) doit être joint.

Sont exclus de la garantie:

- Abrasion habituelle
 - Applications inappropriées, comme par exemple la surcharge de l'équipement, accessoires non certifiés et accessoires modifiés.
 - Dommages causés par des effets étrangers, l'usage de la force ou par des corps étrangers
 - Dommages, qui résultent de la négligence du manuel d'utilisation, par exemple la connexion à une mauvaise tension ou la négligence des instructions de montage
 - Dispositifs complètement ou partiellement démontés
17. L'appareil n'est pas destiné à être utilisé dans un environnement enrichi en oxygène, un environnement d'anesthésie.
18. L'appareil doit être conforme à la norme CEI 60601-1-2 électromagnétique.
19. Ne pas enrouler le câble autour de votre cou pour provoquer l'asphyxie
20. Si l'équipement ne fonctionne pas, vérifiez SVP si la connexion est ferme premièrement, d'autres questions supplémentaires, envoyant au point de garantie est recommandé
21. La modification non autorisée de l'équipement peut entraîner des dangers.

Illustrations (Continued)

Illustration 9.



GlobTek® Inc.
"your power partner"

EN

Operating Instructions

Universal Plug-In Power Supply

Model: GT*41134*****,GT*96060*****,
GTM96060-0706-1.0

The following international GlobTek plug adapters are available:

 US • CAN • JAP NEMA 1-15P / IEC PLUG A (works in Plug B) Configuration: NA 2 Blade, Class II US/CA/JP P/N : Q-NA(R)			

© 2018-2024 GlobTek, Inc. Manufacturer: GlobTek, Inc. 186 Veterans Drive, Northvale NJ 07647, USA. P/N: 809-PSUQKITUSERMANUAL Rev. 4.0. The acceptance of this user manual and use of the product described in this document indicates the acceptance and binding of the customer to GlobTek terms and conditions, which supersede all other agreements, terms, and conditions (<https://en.globtek.com/terms-and-conditions/>). The user of this product (Customer) agrees that they are a sophisticated user who understands the technical specifications and safety implications of use of this product and it is their sole responsibility to determine suitability/fitness of the product described in this specification for their application and to take precautions for protection against any hazards attendant to the handling and use of the product.

Illustrations (Continued)

Illustration 10.

INTENDED USE


The product serves as the power supply for electrical devices. The power consumption by the devices that are connected must not exceed the output current of the power supply. It is imperative that the safety instructions and all other information contained within these operation instructions are adhered to. This product complies with statutory, national and European requirements. All of the company names and product designations are trademarks of the respective owners. All rights reserved.



INSTALLATION


Installation of an INTERCHANGEABLE ADAPTER on the POWER SUPPLY.

- Push the upper edge of the **INTERCHANGEABLE ADAPTER** into the **POWER SUPPLY** at an angle of 30-60° (*image 1*). The upper edge is straight and the lower edge is a horseshoe shape. The power supply is shaped accordingly.
- Now press the adapter into the power supply in a downward movement (*image 2*). The adapter locks in place and makes a clicking sound.
- Now check that the interchangeable adapter has been installed correctly: Hold the power supply in one hand and pull up the adapter with the other hand. The adapter must be fitted securely in the power supply and must not come loose.

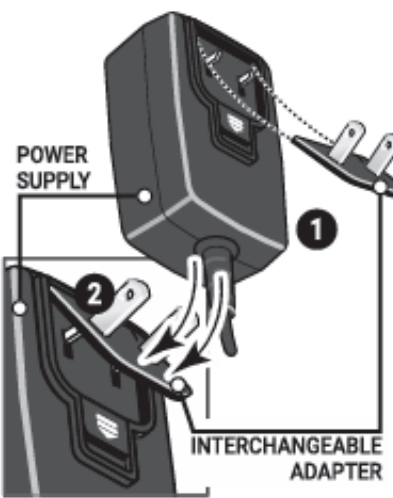

Removing the interchangeable adapter from the power supply

- Using your thumbs or a different finger, push the **FLAP** downward (according to the **ARROW MARKING**) (*image 3*).
- Hold the **FLAP** securely at the bottom and pull the **INTERCHANGEABLE ADAPTER** upward.

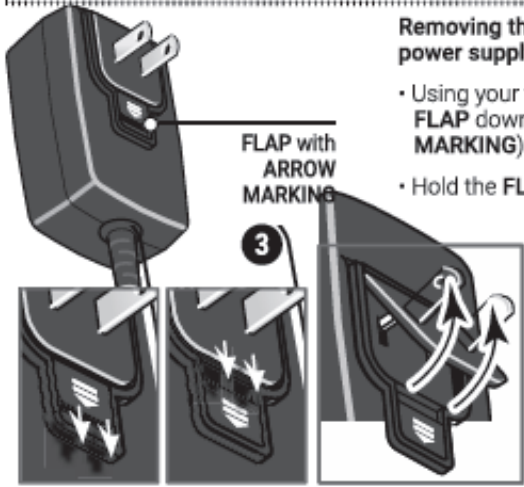

EC

REP

GlobTek Deutschland GmbH
Hafenweg 26a, DE-48155 Münster, Germany
+49 251 13496371 / vertrieb@globtek.de

Illustrations (Continued)

Illustration 11.

SAFETY INSTRUCTIONS



The warranty/guarantee is voided in the event damage is caused by not adhering to the instructions provided in this document. We do not accept any liability for consequential damages. We do not accept any liability for property damage or personal injuries caused by improper handling or not adhering to the following these safety instructions. The warranty/guarantee is voided for the following:

- For safety and approval reasons, (CE) independently converting and/or modifying the product is not permitted . Do not dismantle this product for any reason!
- KEEP AWAY FROM CHILDREN. This device processes electronic current and can be harmful if not handled properly. Keep the product and product packaging away from children.
- KEEP AWAY FROM WATER, MOISTURE AND WATER VAPOR. Use and handle this product in dry and enclosed indoor areas ONLY. DO NOT handle with moist or wet hands.
- This product has been built according to protection Class II and is supplied with international plug adapters. Only use this product in safe and _____-certified power sockets according to regional laws.
- AVOID CONDENSATION. Allow this product enough time to achieve room temperature before installation. Condensation can accumulate with excessive temperature variation and can destroy sensitive electronic equipment and lead to electric shock..
- Avoid the following adverse ambient conditions at the point of installation and during transportation :
 - Dust, flammable gases, vapors and/or solvents
 - Strong impacts, vibrations, and/or shocks
 - Strong magnetic fields that can exist around machines or loudspeakers.
 - Moisture or too high a humidity
 - Coldness or heat, direct sunlight
- Do not pull the electrical cord to disengage product from power source. Hold on either side of the plug housing and pull from the socket.
- Use this product in a well-ventilated area. When operating this product, ensure there is sufficient air circulation and do not cover the product.
- Avoid injury. Always lay the cable out of pathways to avoid tripping injuries.
- Ensure the low voltage cable of the plug-in power supply does not get crushed, twisted or damaged due to sharp edges.
- Avoid fire or life-threatening electric shock. Do not allow product to come in contact with water or liquids. Do not place product within the vicinity of water or liquids (e.g. faucets, vases, fish tanks) to avoid accidental contact.
- Do not touch power supplies or cords if they exhibit signs of damage. Switch off the main power source to the socket from the circuit breaker or fuse box. Once power has been switched off, remove the plug from the socket. Dispose of the damaged device in accordance with environmental regulations and do not use this product further.
- DO NOT leave packaging or other product elements accessible to children.
- Handle with care. This product can be damaged by impacts, shocks, vibrations or falls from short heights.

Illustrations (Continued)

Illustration 12.

EU CONFORMITY



The manufacturer hereby declares that the product conforms to the relevant European directives.

START UP

Important! The output voltage of the power supply must correspond to the input voltage of the electrical device that is connected to the power supply. If this is not adhered to, the power supply and the connected devices will be destroyed. This voids the warranty/guarantee.

Please install the corresponding country adapter that is appropriate for the power socket. When connecting to the device, ensure that the polarity of the plug-in connector is correct. If the device has an on/off switch, this is how you switch on the device. If this does not function correctly, switch it off immediately and check that the cable connection is correct (polarity).

MAINTENANCE & CLEANING

You do not need to perform maintenance on the plug-in power supply. Do not ever dismantle it. Remove the power supply from the power socket prior to cleaning. A dry, soft, clean cloth is sufficient for cleaning. Do not use any aggressive cleaning agents as this can result in discoloration. Do not press too heavily on the surface when cleaning. This helps avoid scratch marks.

DISPOSAL

Dispose of the product at the end of its service life in accordance with the applicable statutory provisions. Within the European Union, the symbol at right indicates this product must not be disposed of in domestic household waste. Old devices contain valuable recyclable material that should be re-used in order to prevent the environment and human health from being damaged due to unregulated waste disposal. For this reason, please dispose of old devices via appropriate collection systems or send the device for disposal at the place where you bought it. This is where the device's materials will be sent for recycling.



STATEMENT of WARRANTY

In accordance with statutory warranty claims, the manufacturer provides a warranty according to the laws of the respective country. The warranty is always for at least 1 year (2 years in Germany). The warranty begins as of the date on which the device was sold to the end consumer. The warranty only covers defects attributable to material or manufacturing errors.


The following cases are excluded from the warranty:

- Normal wear
- Improper uses such as overloading the device or unapproved accessory parts
- Damage due to third party interference, use of force or foreign objects
- Damage arising from non adherence to the instruction manual such as connecting to the wrong mains voltage or non adherence to the assembly instructions.
- Fully or partially dismantled devices

Legal notice: GlobTek Inc., 186 Veterans Drive, Northvale, New Jersey 07647, USA
All rights reserved.

Illustrations (Continued)

Illustration 13.

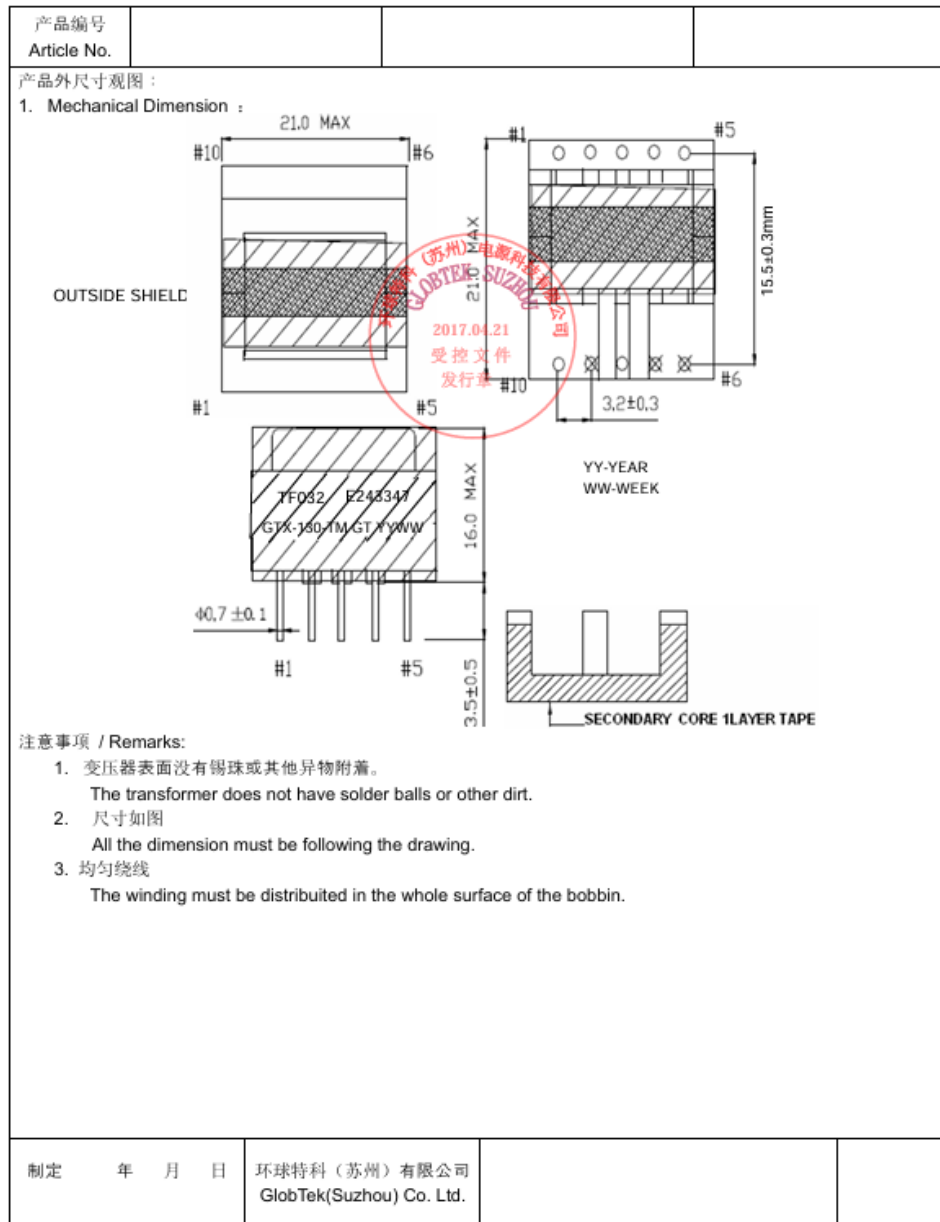

Illustrations (Continued)

Illustration 14.

产品编号 Article No.			
3. 绕线规格: Winding Information:			
顺序 Order	PIN 脚 PIN No.	铜线 Copper wire	圈数 Turns
W1	1---4	2UEW 0.16*4mm Class B	11
2turns W=8.5mm, T=0.025mm Insulation tape			
W2	(3)---5	2UEW 0.13mm Class B	156
2turns W=8.5mm, T=0.025mm Insulation tape			
W3	4---	Copper foil 7mm*0.025mm (非自粘)	1.1T
2turns W=8.5mm, T=0.025mm Insulation tape			
W4	(10)---(8)	TRWB 0.55mm	9
2turns W=8.5mm, T=0.025mm Insulation tape			
W5	2--4	2UEW 0.16mm Class B	11
2turns W=8.5mm, T=0.025mm Insulation tape			
2turns W=4.5mm, T=0.025mm tape Fix the core			
	4--	Copper foil 7mm*0.05mm (非自粘)	1.1T
2turns W=8.5mm, T=0.025mm Insulation tape			
2turns W=14mm, T=0.25mm Insulation tape			
1 2 3 4 5 10 9 8 7 6 注意事项(Note): 1. PIN6,7,8,9 CUT OFF			
PRI SEC W5 W4 W3 W2 W1 GAP NOGAP			
制定	年 月 日	环球特科 (苏州) 有限公司 GlobTek(Suzhou) Co. Ltd.	

Testing Considerations

Samples of the medical power supply were subjected to the following test program with satisfactory results. All tests were conducted in accordance with ANSI/ AAMI ES60601-1:2005, ES60601-1:2005/AMD1:2012, ES60601-1:2005/AMD 2:2021; CAN/CSA C22.2 No. 60601-1:14 + A2:22 (R2022) Medical electrical equipment — Part 1: General requirements for basic safety and essential performance, IEC 60601-1-6 Edition 3.2 2020-07; CAN/CSA-C22.2 NO. 60601-1-6:11 + A1:15 + A2:21 (R2021) Medical electrical equipment - Part 1-6: General requirements for basic safety and essential performance - Collateral standard: Usability, ANSI/ AAMI HA60601-1-11:2015 [Including AMD1: 2021], CSA C22.2 NO. 60601-1-11:15 (R2020) + A1:21 Medical Electrical Equipment -- Part 1-11: General requirements for basic safety and essential performance -- Collateral Standard: Requirements for medical electrical equipment and medical electrical equipment and medical electrical systems used in the home healthcare environment.

Only these tests were considered necessary due to engineering considerations. Detailed test results are on file at MET Laboratories under project number 137720.

Tests Conducted

ANSI/AAMI ES60601-1, CAN/CSA-C22.2 No. 60601-1

- 1 Power Input Test (cl. 4.1.1)
- 2 Humidity Test (cl. 5.7)
- 3 Marking Legibility Test (cl. 7.1.2)
- 4 Marking Durability Test (cl. 7.1.3)
- 5 Leakage Currents And Patient Auxiliary Currents Test (cl. 8.7)
- 6 Dielectric Voltage Withstand Test (cl. 8.8.3)
- 7 Ball Pressure Test (cl. 8.8.4.1)
- 8 Temperature Test (cl. 11.1)
- 9 Overflow, Spillage, Leakage, Cleaning, Sterilization And Disinfection, Harmful Ingress of
- 10 Interruption of Power supply (cl. 11.8)
- 11 Hazardous Conditions And Fault Conditions – Electrical Single Fault Condition (cl. 13.2.2)
- 12 Mechanical strength (cl.15.3)

ANSI/AAMI HA60601-1-11, CAN/CSA C22.2 NO. 60601-1-11

- 1 Environmental conditions of transport and storage between uses
- 2 Continuous operating conditions
- 3 Additional requirements for ingress of water or particulate matter into me equipment and me
- 4 Requirements for mechanical strength for non-transit-operable ME equipment

Conclusion

The product covered by this report have been tested, examined, and found to comply with the applicable requirements of Standards listed on the title page of this report regarding the risks of electrical hazards, mechanical hazards, radiation hazards, excessive temperature hazards only. The risks associated with other properties of the product(s) have not been evaluated. This certification has been granted under a System 3 program as defined in ISO/IEC 17067.

Prepared By:



Jack Gan
Project Engineer,
Partner or Safety Laboratory

Reviewed By:



Jackie Zhao
Project Reviewer,
Safety Laboratory